

Novel Buoyancy Based Cell Selection: X-BACS[™] Technology


Dalip Sethi^{1*}, Jon Ellis², Zelenia Contreras², John Perea², Stephen Truong², and Jillian Miller¹

¹Cesca Therapeutics, Inc., Rancho Cordova, CA; ²ThermoGenesis Corp., Rancho Cordova, CA

Introduction


New developments in cell biology and genetic engineering have revolutionized the field of medicine with the greatest impact observed in the field of Immuno-oncology. Efficient isolation of desired cell populations is the cornerstone in the development of manufacturing and quality control processes for emerging cell therapies. We have developed a buoyancy based cell selection method that isolates desired cell populations from a mixture of cell types, such as mononuclear cell (MNC) fraction preparations from whole blood. Using this method, we have efficiently recovered highly pure populations of T-cells (95%) with good yield (85%). The method is simple to execute with standard laboratory equipment.

Principle of Operation

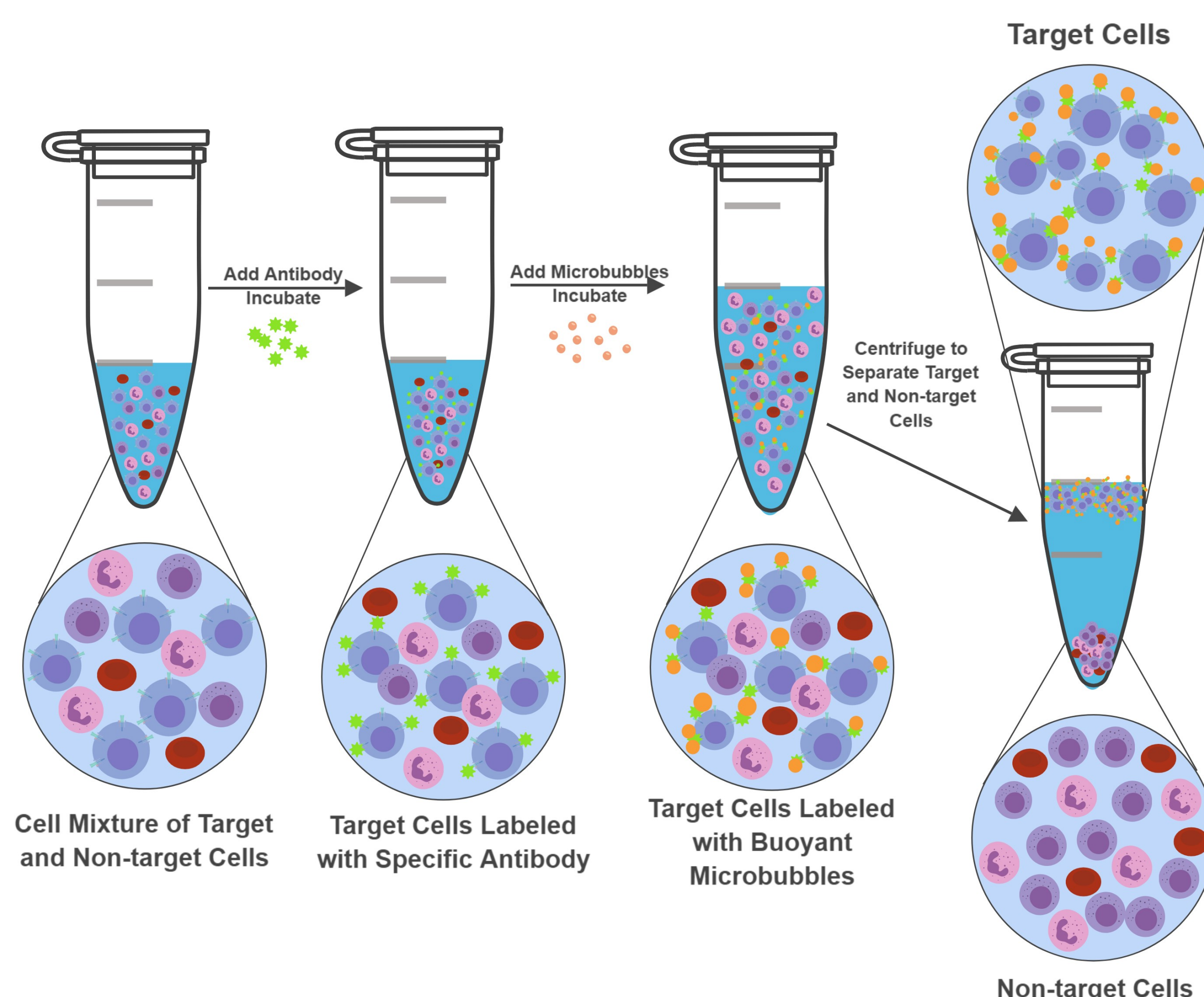

Materials and Methods

Peripheral Blood (PB) units were purchased from BloodSource, Mather, CA. MNC fractions were prepared using both Ficoll-Hypaque density gradient centrifugation or, as an alternative, the X-LAB[™] System.

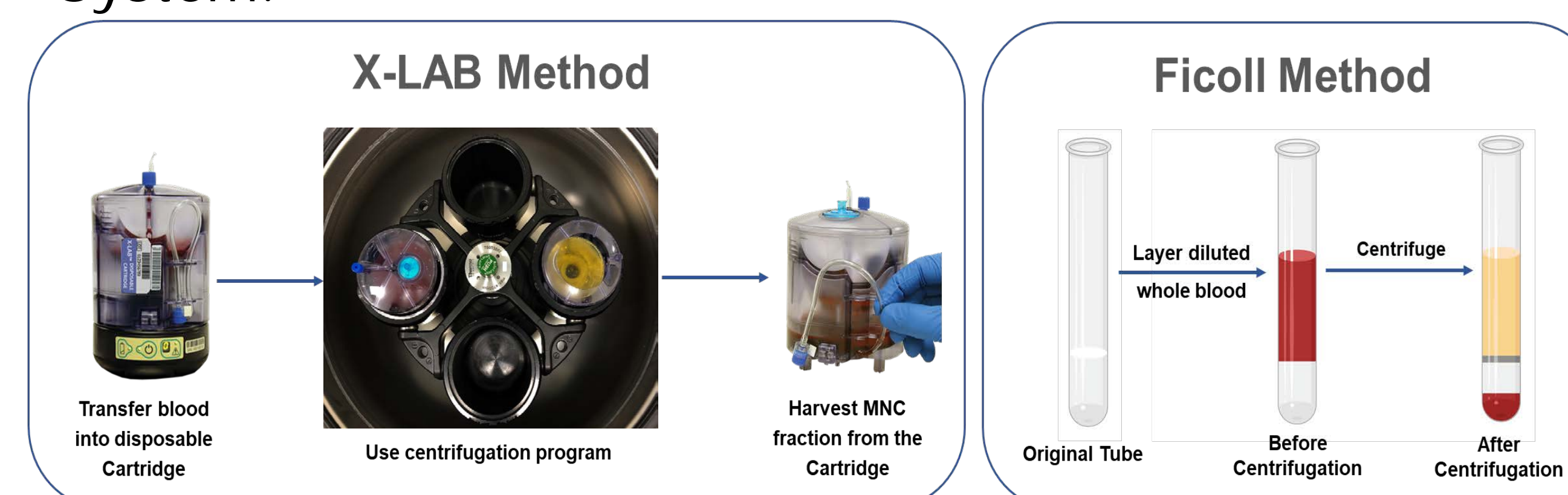

Both MNC fractions were used to select CD3⁺ T-Cells using BACS technology. MNC fractions were incubated with CD3 antibody for 30 minutes, followed by microbubble reagent for 20 minutes. Following incubation, the target and non-target fractions were separated using centrifugation.

Results

Whole blood units (n=4) were split among three (3) independent users and each recovered MNCs using both Ficoll-Hypaque gradient centrifugation as well as with the X-LAB System.

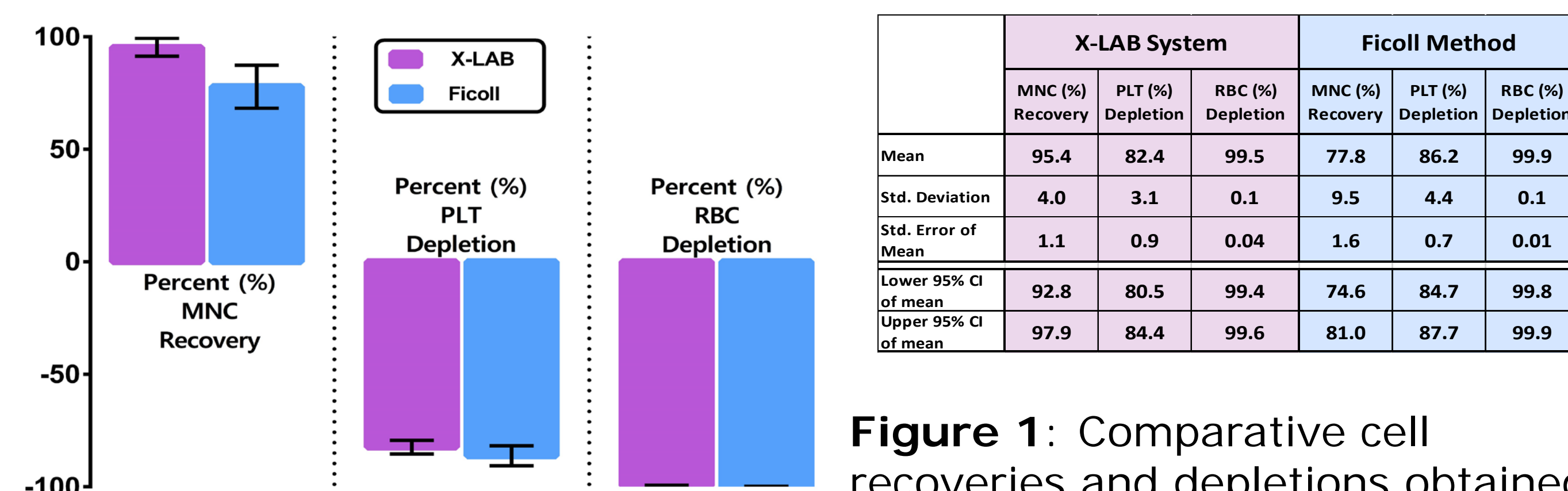


Figure 1: Comparative cell recoveries and depletions obtained using X-LAB Technology and Ficoll method.

CD3⁺ T-cells were selected from both MNC fractions using the X-BACS System. A mean CD3⁺ recovery of 83.5%, 96.3% viability and 95.1% purity was obtained from Ficoll MNC fractions while a mean CD3⁺ recovery of 89.0%, 98.3% viability and 97.3% purity was obtained from X-LAB MNC fractions.

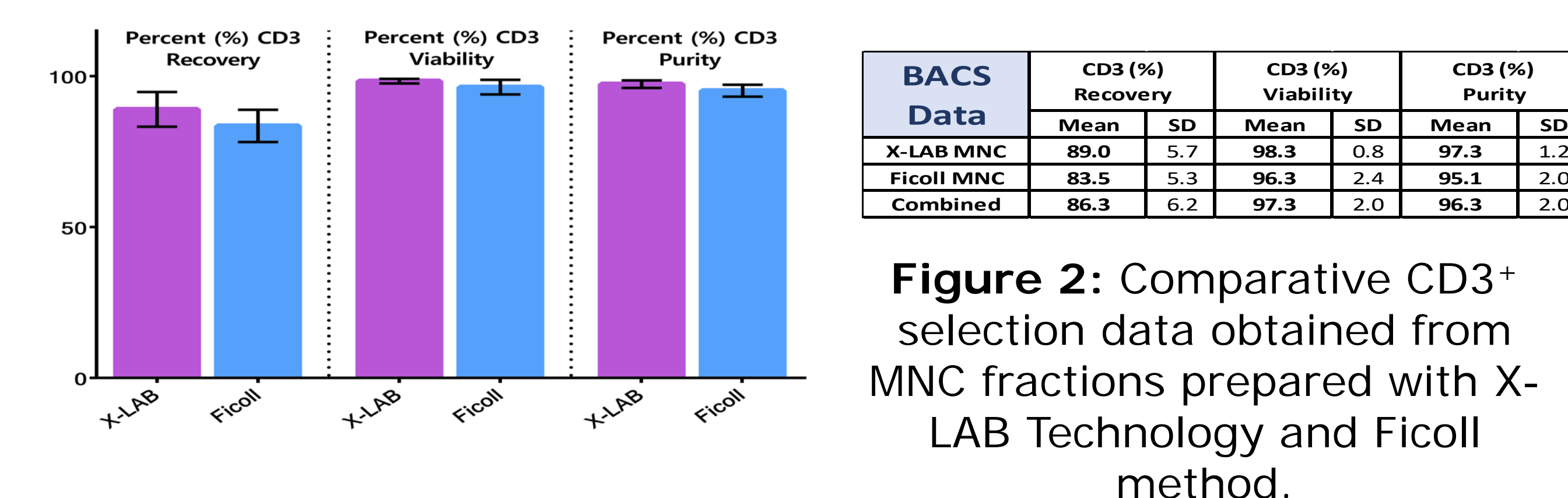


Figure 2: Comparative CD3⁺ selection data obtained from MNC fractions prepared with X-LAB Technology and Ficoll method.

The overall percent (%) CD3⁺ cell recovery was calculated by multiplying percent (%) MNC recovery with percent (%) CD3⁺ BACS recovery. In summary, better CD3⁺ selection was obtained using the X-LAB System followed by BACS compared to Ficoll method followed by BACS.

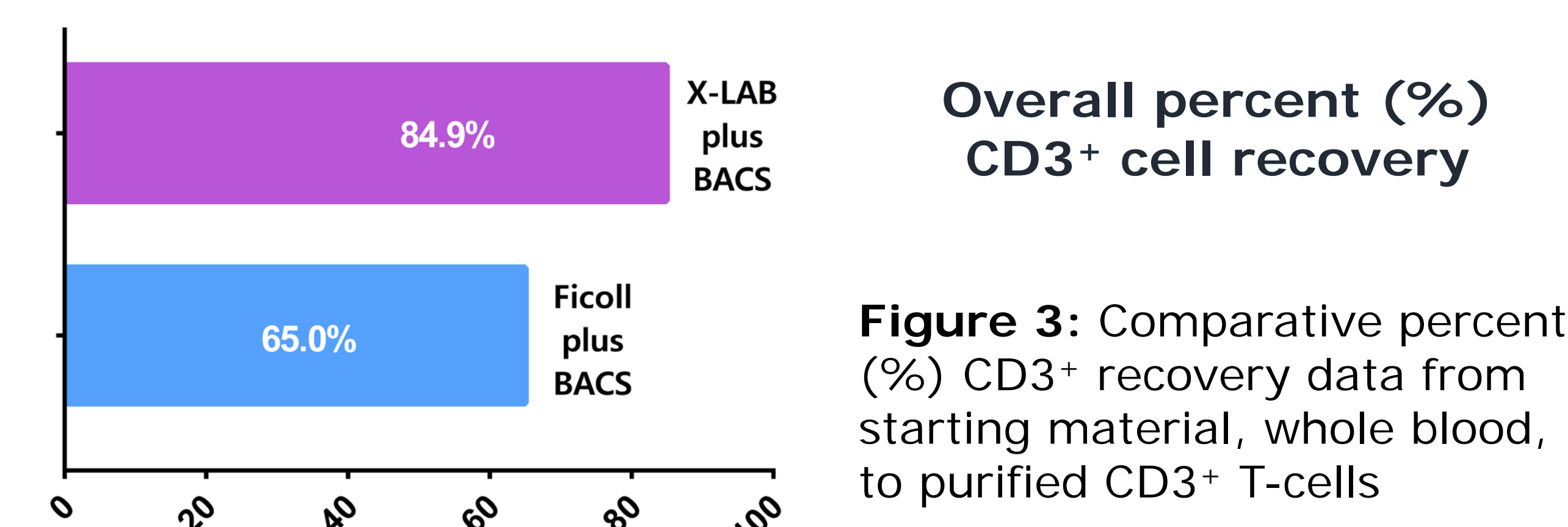


Figure 3: Comparative percent (%) CD3⁺ recovery data from starting material, whole blood, to purified CD3⁺ T-cells

Conclusions

We have developed a buoyancy based cell selection method that provides an efficient means of isolating a highly pure population of CD3⁺ T-cells (>95%) with good yield (>85%) and excellent viability (>95%).