




X-WASH® System

The X-WASH System washes, resuspends and concentrates cell suspensions within a sterile single-use disposable.


- 
- Single-Use, Sterile, Functionally Closed System
 - Reliable High Recovery (>80%) and Viability (>90%) of Cells
 - Process 240 mL down to 3mL in 10-15 minutes
 - Centrifugation-based Cell Washing System
 - GMP-compliant System
 - Customized Protocols

X-WASH® System

Semi-automated Cell Washing with the X-WASH System

For ordering information, contact your local distributor or visit us at www.ThermoGenesis.com


- 1** Load cells into Disposable Cartridge

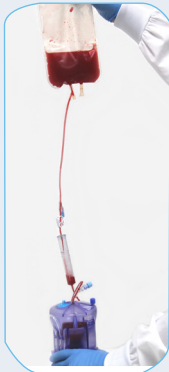

- 2** Place Cartridge and Control Module in Centrifuge



- 3** Retrieve Washed Cells

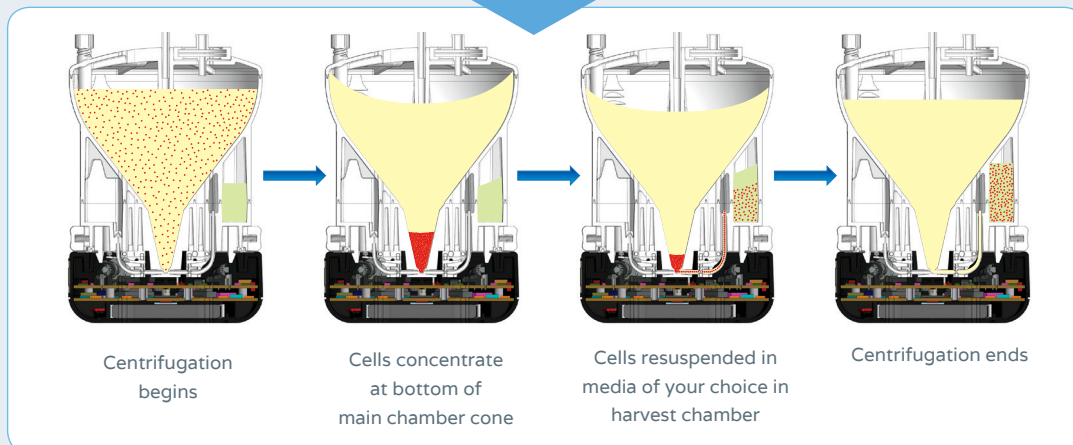

The ThermoGenesis Advantage

Feature	X-WASH System
Customized Protocols	Yes
Process Volume	40 to 240 mL
Harvest Volume	Programmable (3 – 40 mL)
Average Process Time	10-15 minutes/wash cycle
GMP Compliance	Electronic transfer of process data and barcode compatible
Cell Recoveries	>80%
Closed System	Yes

ThermoGenesis Corp.

2711 Citrus Road, Rancho Cordova, CA 95742

Toll-free: 800.783.8357 (US/Canada)

Direct: 916.858.5100

Fax: 916.858.5199

